

Data Sheet

Anti-FLAG Tag

Mouse Monoclonal Antibody
Anti-FLAG Tag (FG4R)
Catalog #: w35014

Specificity:

Recognizes the N-terminal, C-terminal or internal FLAG-tagged fusion proteins

Immunogen:

DYKDDDDK (FLAG) synthetic peptide conjugated to KLH

Applications:

Dot Blot, ELISA, IP, ICC, WB

Image:

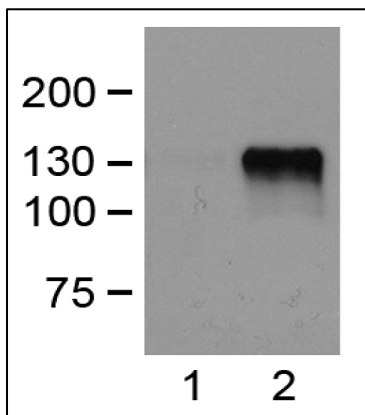


Fig. 1:1000 (1 μ g/mL) Ab dilution probed against HEK293 cells transfected with FLAG-tagged protein vector; untransfected (1) and transfected (2).

Clonality | Host Species | Isotype:

FG4R | Mouse | IgG2b

Formulation:

1mg/mL in 10mM PBS, pH 7.2, with 0.05% NaN₃

Sensitivity:

WB (with ECL): 1-5ng of FLAG-tagged fusion proteins

Working Conditions:

WB (with ECL): 1:1000-3000 dilution (incubate for one hour at room temperature)
ICC: 1:500-2000 dilution
For best results with other assays (e.g.: Dot Blot, ELISA, IP, etc.), please determine optimal working dilution by titration test

Purification:

Protein A affinity chromatography from mouse ascites fluid

Storage:

Store at -20°C for up to 2 years. Centrifuge after first thaw to maximize product recovery. Aliquot to avoid repeated freeze-thaw cycles. Store aliquots at 4°C for 1-2 weeks

OUR PRODUCTS ARE FOR IN VITRO RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.