

Data Sheet

Caspase-7

Human, recombinant, C-terminal His Tag
Catalog #: w80011
Lot#: 130116 **Conc.:** 2.28 mg/ml

Formulated in: 40 mM Tris-HCl, pH 8.0,
110 mM NaCl, 2.2 mM KCl, 155 mM
imidazole, 20% glycerol and 3 mM DTT.

Stability: >6 months at -80°C

Reference(s):

1. Wang DH., et al., PLoS one, **7(6)**:e39920.
doi:10.1371/journal.0039920. (2012).


Description: Human Caspase-7 (full length) (GenBank Accession No. NM_001227) with C-terminal His tag, MW= 35.1 kDa, expressed in an E. coli expression system. Active caspase 7 is a heterotetramer consisting of two copies each of ~20 kDa and ~11 kDa subunits.

Specific Activity: 5,337 pmol/min/μg.
Assay conditions: HEPES, pH 7.4, 10 mM EDTA, 0.05% Chaps, 5 mM DTT, and 10 μM Caspase substrate. Incubate for 30 min at RT. Fluorescence intensity measured at exc. 380 nm, em. 500 nm.

Application: Useful for the study of enzyme kinetics, screening inhibitors, and selectivity profiling.

Quality Assurance

Specific Activity

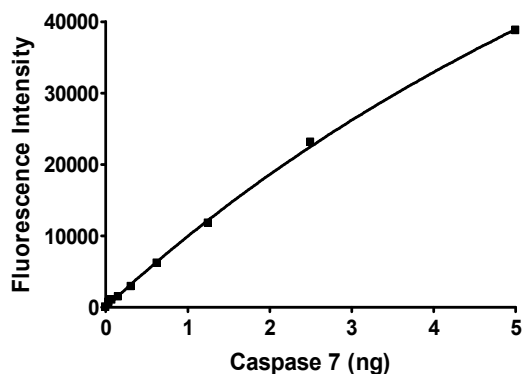


4-20% SDS-PAGE
Coomassie staining

Lane 1:
4.56 μg Caspase-7
(subunits 1 & 2)

Lane 2:
Protein Marker

MW: 11 kDa
and 20 kDa

Purity: ≥95%

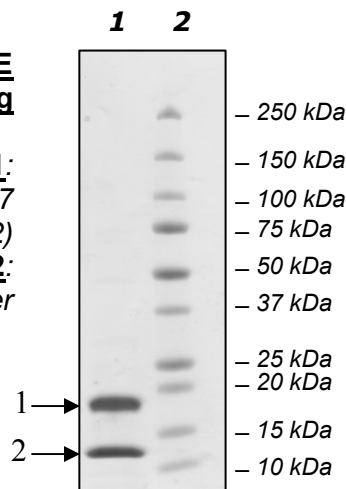

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.



9891 Irvine Center Dr. Suite 200
Irvine, CA 92618
United States
Tel: 1.800.831.1518
Fax: 1.800.831.1518
Email: sale@westbioscience.com

Protein Sequence:

MADDQGCIEEQGVEDSANEDSVDAKPDRSSFVPSLFSKKKKNVTMRSIKTTTRDRVPTY
QYNMNFELGKCIINNKNFDKVTGMGVRNGTDKDAEALFKCFRSLGFDVIVYNDSCA
KMQDLLKKASEEDHTNAACFACILLSHGEENVYIGKDGVTPIKDLTAHFRGDRCKTLLE
KPKLFFIQACRGTELDGIGADSGPINDTDANPRYKIPVEADFLFAYSTVPGYYSWRSP
GRGSWFVQALCSILEEHGKDLEIMQILTRVNDRVARHFESQSDDPHFHEKKQIPCVVS
MLTKELYFSQH

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.