

Data Sheet

Cathepsin B

Human, recombinant, C-terminal His-tag

Catalog #: w90012

Lot #: 111101

Conc.: 2 mg/ml

Formulated in: 40 mM Tris-HCl, pH 8.0,
110 mM NaCl, 2.2 mM KCl, 250 mM
imidazole, 20% Glycerol

Stability: >6 months at -80°C

References:

1. Tu, C. *et al.*, *Cancer Res.* 2008 Nov 15;**68(22)**:9147-56.
2. Hwang, JH. *et al.*, *Cancer Lett.* 2009 Jan 18;**273(2)**:266-72.

Description:

Human cathepsin B (GenBank Accession No. NM_001908), CD33 signal peptide (1-16 a.a.) + Cathepsin B (18-339 a.a.), with C-terminal His-tag, MW=43 kDa, expressed in FreeStyle 293-F cells.


Specific Activity: 4048 pmol/min/μg

Assay condition: 25 mM MES, pH 5.0, 10 μM substrate (Z-Leu-Arg-AMC), 10 min, room temperature. Prior to assay enzyme was activated with preincubation at 10 μg/ml in 25 mM MES, pH 5.0, 5 mM DTT for 15 min at RT.

Application:

Useful for the study of enzyme kinetics, screening inhibitors, and selectivity profiling.


Quality Assurance

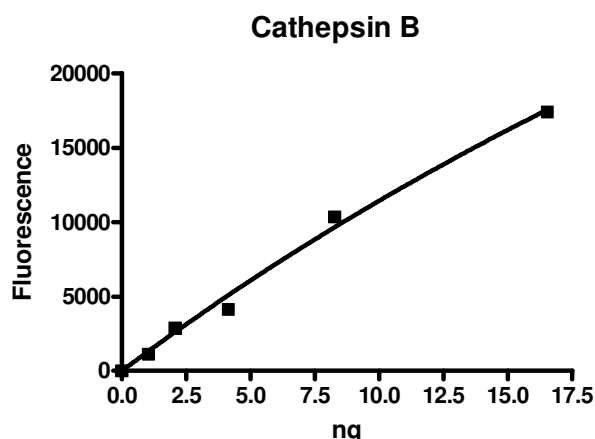

10 % SDS-PAGE Coomassie staining

Lane 1:
7 μg Cathepsin B
Lane 2:
Protein Marker
BioLabs (#P7708L)

MW: 43 kDa
Purity: ≥95%

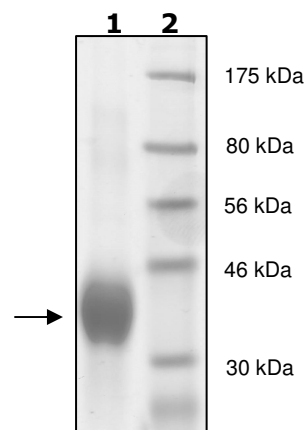

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.